ES

WKB006  
Manual de instrucciones  
WE-2020-WKB006-1

# TERMOTERAPIA

## ADVERTENCIA

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# TERMOTHERAPY

## WARNING

- Never under any circumstances attempt to open or inspect the internal components or accessories of the machine. If an inspection becomes necessary, please contact your supplier or a suitably competent technician.
- Never use the unit in humid or poorly ventilated environments, such as wet rooms, bathrooms, etc. to prevent damage to the electrical components of your unit.
- Never use the unit outdoors; exposure to the elements can cause damage.
- If the power cord was destroyed, please ask a professional to service and replace it to avoid hazards, or contact the dealer for service.

## EC-DECLARATION OF COMPLIANCE

This product meets the following requirements:

1. Conditions set out in the Low Voltage Directive: 2014/35/EC
2. Conditions set out in the Electromagnetic Compatibility Directive: 2014/30/EC

## RECYCLING

Whenever any of our electrical products or their presentation packaging bear this symbol, it means that these products should not be treated as conventional household waste in Europe. To ensure proper treatment of this waste, dispose of it in accordance with local laws or as required for the disposal of electrical equipment. This initiative will help preserve save natural resources and

improve standards of environmental protection in regard to the treatment of electrical waste.

## WELCOME!




Thank you very much for purchase. Please read the instructions carefully before operating the instrument. You can also ask a professional for guidance to make sure you are using the instrument correctly. We will accept no responsibility for any accidents that occur due to incorrect operation.

The company reserves the right to amend the details of the instruments with no prior notice. If any mistakes are found in these instructions, please feel free to send us an e-mail informing us so that we can correct them.

## PACKAGE CONTENTS

- These components are subject to change without prior notice.
- Appearance of product may differ from the image.

Nr.	PRODUCT	REF.	QTY.
1	THERMAL BLANKET	40002.1	1
2	MAIN UNIT	-	1
3	FEEDER CABLE	-	1

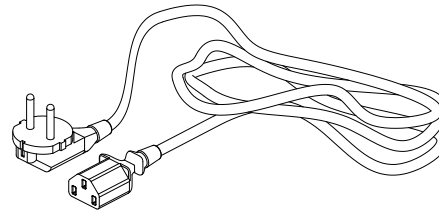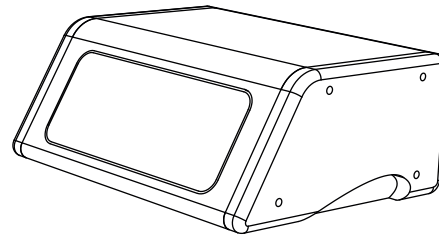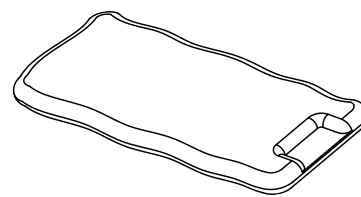4. This device has three independent and adjustable temperature zones.

5. First adjust the treatment time, from 1 to 60 minutes, by using the buttons +/-.

6. Adjust the operating temperature for every zone by using the buttons +/- and press "START / PAUSE" in each channel. The adjustment range is from 40 to 85°C. (The thermal comfort range is from 30 to 50°C.)

7. To turn off the device, press again the "START" button on the right.

## CAUTIONS

- If the surface of the blanket is broken or aging, please replace it with a new blanket to prevent accidental skin burns or electric shocks.
- The heating blanket can generate light magnetic radiation to the human body due to its high power. Please note that pregnant women, children, elderly people and people with sensitive skin should avoid using this instrument. Please drink plenty of water before treatment to avoid dehydration and heatstroke.

- Please check with the patient if the operating temperature is appropriate during treatment.


## MAINTENANCE

- The leather surface of this blanket is made of PVC material. Please keep it dry and clean after use. If the leather is stained with weight-loss cream or other oily substances, use a neutral cleaning liquid to wipe them off. Never use acid or alkaline organic liquids for cleaning.
- Please cut off the power if the instrument will not be used for a long time. Store the housing and accessories separately in a dry environment.


## ACCESSORIES

REF.	DESCRIPTION
40002.1	THERMAL BLANKET

## FRONT PANEL

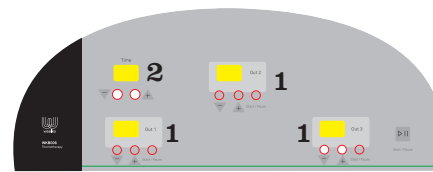1. Temperature control for the three zones
2. Timer

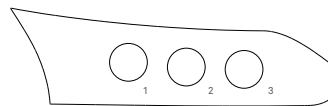3. Air tubes output

## INSTRUCTIONS OF USE

We recommend you to give a massage before and after the thermotherapy treatment.

1. Connect the power cables to the rear panel, as indicated in the sketch as "POWER LINE INPUT".
2. Connect the 3 accessories properly to the front panel, as indicated in the sketch.
3. Press "START" button on the right.

## TECHNICAL SHEET

Weights and sizes	Packing size (cm)	58x41x52
	Gross weight (kg)	9
	Product size (cm)	48x13x32
	Blanket size (cm)	180x80
Working characteristics	Working power	300 W
	Temperature range	30 - 50°C
	Maximum temperature	55°C
	Preset working time	30 min
Electrical characteristics	Rated voltage	220V ~ 240V
	Frequency	50Hz/60Hz

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